


ACHIEVE HIGH SURVIVAL RATES WITH The Cryotec Method

Cryotec products are designed for safe and effective operability

Plates
designed for
safer embryos
manipulation



Cryotec
optimized for
ease of
operation

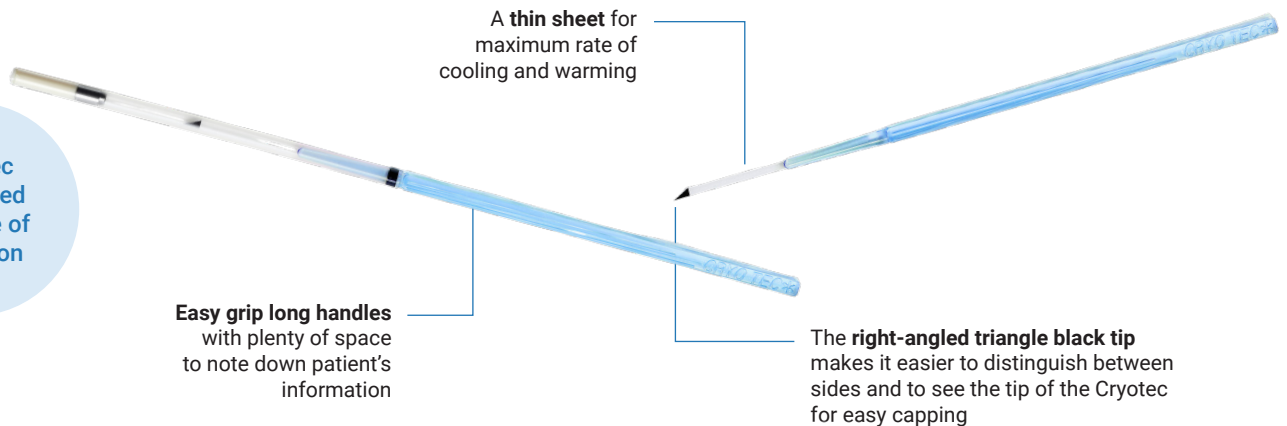

Solutions enable safe recovery

- Safe and long-lasting due to solutions being free from HSA and proteins
- Low impact on cells due to the use of endotoxin-free trehalose and a low concentration of CPA
- Hydroxypropyl cellulose (HPC) increases viscosity and slows down the movement of embryos



The Cryotec Method is a **safe and effective** technique for vitrification and warming of oocytes and embryos. It offers a unified protocol for oocytes and embryos of all stages leading to **high survival rates**.



**NO TIME
PRESSURE
DURING THE
PROCEDURE**



**GENTLE
OSMOTIC
CHANGES**



**SAFE
SOLUTIONS
WITH RELATIVELY
LOW TOXICITY**



**EFFECTIVE FOR
SENSITIVE
OOCYTES/
EMBRYOS**



**KEEPING
FERTILITY
OPTIONS
OPEN**

Our Products

Vitrification Kit 101 (For 3 applications)



1.0 ml Equilibration Solution (ES) x1
1.0 ml Vitrification Solution (VS) x2
Cryotec x4
Vitri Plates x3

Vitrification Solution Set 110 (For 10 applications)



1.8 ml Equilibration Solution (ES) x2
1.8 ml Vitrification Solution (VS) x4

Vitri Plate (Set of 10)



Warming Kit 102 (For 1 application)



1.8 ml Thawing Solution (TS) x1
0.5 ml Diluent Solution (DS) x1
1.0 ml Washing Solution (WS) x1
Warm Plate x1

Warming Solutions Set 205 (For 5 applications)



1.8 ml Thawing Solution (TS) x5
1.8 ml Diluent Solution (DS) x1
1.8 ml Washing Solution (WS) x2

Warm Plate (Set of 10)



Cryotec

10 Cryotec vitrification devices



Available in blue, green, yellow, white and red

How to order

Please contact us by email
globalsales@reprolife.jp or
scan QR code to visit our
contact page



Note: Cryotec products are designed for simple, safe and effective operability to give the best chance for successful vitrification and warming of oocytes and embryos. REPROLIFE Inc. is ISO 13485 certified. We only ship products that have passed a variety of rigorous tests.

Disclaimer: Some products are not available in the USA, EU and other countries. Contact us to inquire.

Produced by REPROLIFE Inc.



Version code: This brochure was
E_102 updated in December 2024

